

N 1500

# ELECTROPORATION

HIGH-POWERED

UNIQUE WITH EXPONENTIAL AND DECAMETRIC PULSED CURRENTS.



**N1500**  
PROFESSIONAL EQUIPMENT

- BIPHASIC PULSE.
- HIGH VOLTAGE AND LOW INTENSITY.
- TOUCH SCREEN DISPLAY.
- IMPEDANCE DETECTOR (ID®).
- BODY AND FACIAL DEVICE.
- EXPONENTIAL WAVES.
- DECAMETRIC WAVES.
- RUNNING-TIME COUNTER.
- SMT MICROELECTRONICS.
- METAL CABINET.

## THE MARKET'S NOVEL DEVICE

In an unstoppable search to satisfy the most demanding professional needs Noto® created N 1500®, combining technology, design, and performance for electroporation treatments. N 1500 includes "ID®" (Impedance Detector) system - DEVELOPED BY NOTO®.

### "ID®" SYSTEM (IMPEDANCE DETECTOR)

ID® smart system reveals the professional, with a luminous sensor, the accurate patient-electrode link and the correct emission of the selected current. This new technology assures treatment safety to professionals.

### MICROELECTRONICS TECHNOLOGY

N 1500 is developed under strict safety and quality standards, like each Noto device and development. Its hardware and manufacturing system (with SMT technology) gives N 1500 excellent usage, stability, and reliability features for the patient and professionals while treatment is in progress.



NOTO proposes a different course for aesthetics, a new scenario where therapies combination, cutting-edge technology, and quality products and supplies are fundamental for visible, effective, and durable results - with the launched of N product Line.

## ELECTROPORATION

Electroporation is the non-invasive alternative to conventional mesotherapy.

Applies short electrical pulses acting on cells membranes increasing its permeability to facilitate active ingredients passage. This technique is up to 350 times more effective than ordinary application.

N 1800 allows transdermic transference of micro or macromolecules greater than 800.000 Dalton (Da), as hyaluronic acid or heparin - or 500.000 Dalton as collagen.

### MAIN ADVANTAGES

- Biphasic device.
- Painless and doesn't produce secondary effects.
- Friendly and easy of use.
- Place drug or active ingredient on the application zone, placing and moving the header on the selected area until its absorption.
- It's cold & heat-less. Athermic current.

### MOST FREQUENT APPLICATIONS


- Cellulite
- Skin Stains
- Wrinkles
- Stretch marks
- Flaccidity
- Skin aging
- Among others.

### TECHNICAL DATA N 1500:

MANUFACTURED BY:	NOTO GROUP S.A.
COUNTRY:	ARGENTINA
NAME OF THE PRODUCT:	N 1500
WIDE:	185 MM
LENGTH:	260 MM
HEIGHT:	130 MM-75 MM
WEIGHT:	2,1 KG
POWER SUPPLY VOLTAGE:	98 to 240 V- 50/60HZ
POWER:	22 VA
PROGRAMMING:	MICROCONTROLLER
DISPLAY:	128X64 GRAPHIC TOUCH
SOURCE:	INTERNAL SWITCH
TECHNOLOGY:	SMT
TIME CONTROL:	MANUAL

ELECTROPORATION SET:	
1 ACETAL RESIN APPLICATOR	
1 BODY SURGICAL STAINLESS STEEL HEADER OF 32 MM DIAMETER.	
1 FACIAL SURGICAL STAINLESS STEEL HEADER OF 15 MM DIAMETER.	
ELECTROPORATION:	
BIPHASIC WAVE FORM:	OWN
FREQUENCY:	2,1 KHZ
OUTPUT VOLTAGE:	800 VPP with 47 K charge
OUTPUT SYSTEM:	BIPOLAR
APPLICATOR HEADER:	SURGICAL STAINLESS STEEL
FACIAL:	15 MM
BODY:	32 MM
APPLICATOR MATERIAL:	ACETAL RESIN
CABLE:	INJECTED WITH WIDE CONNECTOR INPUT

NOTE: DEVICE SPECIFICATIONS MAY CHANGE WITHOUT PREVIOUS NOTICE.

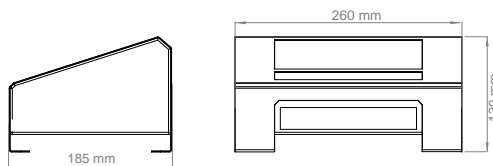

NOTO supports all clients with a set of tools available at all times.

Tools as training, post-sale service support, and assured spare parts.  
All these instruments with the aim of making Noto's products your best choice.

More equipments by NOTO:

N 1400 - MICRODERMABRASION.

N 1104 - BIPHASIC ELECTROTHERAPY. T.E.N.S. + RUSSIAN WAVES + SQUARE WAVES. 26 PROGRAMS.