

NEW. THE BEST COMBINATION

TRICAV

A NEW CONTOURING STANDARD.
ULTRACAVITATION + TRIPOLAR RADIOFREQUENCY.



NOTO... ONGOING INNOVATION.

NOTO® gathered TRIMAX® and CAV47® main features to create a brand-new equipment for benefits, results, and safety on treatments...the TRICAV®.

TRICAV MAXIMIZES CONTOURING TECHNIQUE.

Combination of both techniques enables localized adiposity elimination and a renewed skin from the first session.

Diminish volume, improve cellulite and wrinkles, and stimulate collagen and elastine fibers.

Thanks to its latest Tripolar Radiofrequency + Cavitation, TRICAV® enhances the modeling technique and improves results by nourishing in each session tightening and contouring effects.

TRICAV® sets a new path for skin care and body outline, thanks to its effectiveness, reliability, and results levels.


WELCOME TO NOTO! WELCOME TO TECHNOLOGY!

NOTO^F
TRUSTED TECHNOLOGY

NOTO proposes a different course for aesthetics, a new scenario where therapies combination, cutting-edge technology, and quality products and supplies are fundamental for visible, effective, and durable results.

HOW THIS NOVEL SYSTEM WORKS?

Thanks to emitters quantity, positioning, and its cross fields system made of surgical stainless steel, TRICAV® increases its focused energy on the application zone. Temperatures reaches up to 42°C (107.6°F) on epidermis layer and between 55/58°C (131/136.4°F) in dermis layer. This concentration of energy allows an immediate contraction of collagen and elastine fibers stimulating its unnatural process and achieving instant, visible and concrete results.

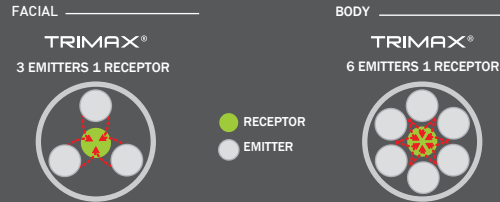

TOP OF THE LINE IN SCULPTING

Ultra-cavitation is a body modeling technique based on micro-bubbles generation. Pressure changes, compression and expansion of the adipocyte's cellular membrane, creating an implosion and the subsequent adipocyte's natural elimination.

TRICAV® creates stability and precision due to its double barium titanate crystal, granting maximum physical contact power of 47 Watts. Work in a friendly and easy way on soft, hard, and edematous subcutaneous fat deposits, with this new equipment.

TRICAV® is effective, with visible results in few sessions, painless, comfortable for the patient, and the best option to fight cellulite.

"ID" SYSTEM (IMPEDANCE DETECTOR)

ID® smart system reveals the professional, with a luminous sensor, the accurate patient-electrode link and the correct emission of the selected current - developed by Noto®. This new technology assures treatment safety to professionals.




TRICAV IPX 0 Device

MANUFACTURED BY:	NOTO GROUP S.A.
COUNTRY:	ARGENTINA
NAME OF THE PRODUCT:	TRICAV
WEIGHT:	4.6 KG.
POWER SUPPLY VOLTAGE:	98 to 240 VAC 50/60HZ
ENERGY CONSUMPTION: RF	120 VA
POWER:	Body: 60 Emission Watt Facial: 40 Emission Watt
CAV POWER:	SHOCK WAVE OF 47 WATT
PROGRAMMING:	MICROCONTROLLER
RF FREQUENCY:	1 MHZ +/- 10%
CAV FREQUENCY:	40 KHZ +/- 10%
DISPLAY:	128X64 GRAPHIC TOUCH
SOURCE:	INTERNAL SWITCH
TECHNOLOGY:	SMT

TIME CONTROL:	MANUAL
ENTRANCE FUSE TYPE:	F 1.5 AL 250 V
EMITTER ENTRANCE FUSE:	F 1 AL 250 V

NOTE: DEVICE SPECIFICATIONS MAY CHANGE WITHOUT PREVIOUS NOTICE.

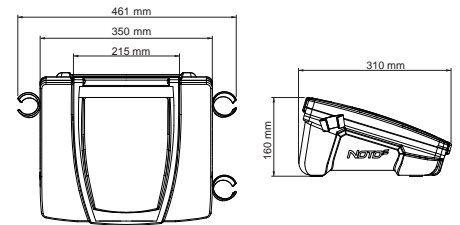

TRICAV

NOTO supports all clients with a set of tools available at all times.

Tools as training, post-sale service support, and assured spare parts.
All these instruments with the aim of making Noto's products your best choice.

NOTO N LINE equipments:

- N 1400** - Microdermabrasion
- N 1500** - High-powered Electroporation