



**N 1104 PROGRESSIVE**

# ELECTROTHERAPY

**RUSSIANS + SQUARE + T.E.N.S. WAVES**

¡100% DIGITAL! ¡54 PROGRAMS!



**2** YEARS WARRANTY  
+ TRAINING

**LATEST TECHNOLOGY**  
PROFESSIONAL DEVICE

- BIPHASIC DEVICE
- 4 INDEPENDENT CHANNELS
- 1 POWER CONTROL PER CHANNEL
- DIGITAL DISPLAY
- START AND END ON RAMP
- PROGRAM CUSTOMIZATION

### THE MARKET'S MODERN EQUIPMENT

N 1104 PROGRESSIVE is 100% digital, biphasic, and with 4 full independent output channels. The device has an isolated transformer in each channel and the particular characteristic of starting and ending treatments on ramp. All these features to achieve safer and reliable devices.

Beyond changing settings at the desired moment, N 1104 PROGRESSIVE has 34 preset programs and 20 customizable programs to create your own parameters and store them. It is presented on a metal cabinet, with a 20 x 4 display.

### MICROELECTRONICS TECHNOLOGY

Its technology, like all developed by NOTO® is manufactured with microelectronics. This characteristic grants excellent stability and reliability usage, enhancing safety, durability, and comfort for the patient.

### A DIFFERENT DEVICE

Manufactured on metal cabinet, entirely digital, silicone injected cables and their identification by colors for quick lecture and application. Different, delivers visible results, quality in treatments, patient comfort, and usage safety are some of N 1104® PROGRESSIVE virtues.

WELCOME TO TECHNOLOGY! WELCOME TO NOTO!

**NOTO**  
THE DIFFERENCE OF BEING UNIQUE

NOTO proposes a different course for aesthetics, a new scenario where therapies combination, cutting-edge technology, and quality products and supplies are fundamental for visible, effective, and durable results.

#### RUSSIAN WAVES

Are modulated medium frequency currents, delivered in low frequency signals. Its medium frequency objective is to decrease skin's impedance to allow low frequency work on muscles. These currents main target is muscle tone or body outline contouring. N 1108 is available to work with bipolar or tetrapolar interferential wave mode.

#### MAIN TREATMENTS:

- Flaccidity and body contouring.
- Atrophy due to inactivity.
- Progressive muscle hypertrophy.
- Glutei hypertrophy and tone.
- Post surgery strenghten process.
- Sports rehab.

#### T.E.N.S.

T.E.N.S. programs utilizes nerve system mechanisms to mitigate chronic and acute pain. Decrease of pain, creating an immediate relief without side effects, is T.E.N.S. main advantage.

#### MAIN TREATMENTS:

- Nerve compression and reflex sympathetic dystrophy (carpal tunnel syndrome).
- Stump and/or phantom member pain.
- Back and neck pain.
- Peripheral Nerves injuries.
- Posttherpetic neuralgia.

#### SQUARE WAVES.

They are presented in low and medium frequency pulses. These waves are the most energetic waves of all electrostimulation currents.

#### MAIN TREATMENTS:

- Body Contouring.
- Sports.
- Tone.
- Muscle firmnes.
- Strength & endurance.
- Contouring.
- Strenghtening.

#### THE DEVICE N 1104 PROGRESSIVE COMES WITH:



4 LARGE Electrodes

4 SMALL Electrodes

6 Elastic Straps

4 Connection Cables

1 User's Manual


#### TECHNICAL DATA N 1104 PROGRESSIVE:

MANUFACTURED BY:	NOTO GROUP S.A.
COUNTRY:	ARGENTINA
NAME OF THE PRODUCT:	N 1104
WIDE:	185 MM
LENGHT:	260 MM
HEIGHT:	130 MM 75 - MM
WEIGHT:	2,1 KG
POWER SUPPLY VOLTAGE:	98 to 240 VAC 50/60 HZ
POWER:	21 WATT
PROGRAMMING:	MICROCONTROLLER
WAVE TYPES:	RUSSIAN, T.E.N.S. and SQUARE 90 MA . 1K
MAXIMUM CURRENTS:	IMPEDANCE
QUANTITY OF CHANNELS:	4 for BODY, ISOLATED INDEPENDENTLY

PROGRAMS:	34 PRESET
TIMER:	1 to 30 MINUTES
CONTRACTION:	2 to 15 SECONDS
RELAXATION:	2 to 50 SECONDS
FREQUENCY:	2 to 150 HZ
RELATION:	1:1 - 2:1 - 4:1
ASCENT RAMP:	1 Second
DESCENT RAMP:	1 Second
SEQUENCE:	1,2,3,4

NOTE: DEVICE SPECIFICATIONS MAY CHANGE WITHOUT PREVIOUS NOTICE.

2500 HZ. RUSSIAN WAVES FREQUENCY CARRIER

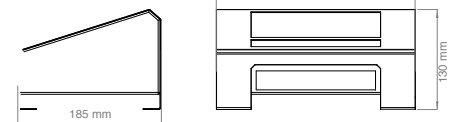

NOTO supports all clients with a set of tools available at all times.

Tools as training, post-sale service support, and assured spare parts. All these instruments with the aim of making Noto's products your best choice.

OTROS EQUIPOS DE LA LINEA N DE NOTO:

N 1800 - MICRODERMABRASION + ELECTROPORATION WITH ID SYSTEM  
CAV 47 - VARIABLE FREQUENCY CAVITATION

